

Specification

Article	Specification
Magnetic Power	3 Tesla
Channel	2 Ch
Magnetic Frequency	1~100Hz
Treatment Protocol	Auto Mode : 20 Manual Mode : 20 User Mode : 20
Treatment Time	1~60min
Interface	8" LCD, Touch Screen, Jog Shuttle
Dimension	541.6 x 500.7 x 1072.9mm
Weight	60kg

The Birth of Salus Talent Pro inspired by 10 years of experience

The World's 1st *Development 3 Tesla,
2Channel Magnetic Field Therapy*

Salus Talent Pro



S a l u s - T a l e n t P r o



Seoul Branch 4F, 224, Baekjegobun-ro, Songpa-gu, Korea
Tel : +82-1588-7395 Fax : +82-2-418-0986

Factory #301~303, 187, Techno 2ro, Yuseong-gu, Daejeon-si, Korea
Tel : +82-42-934-5560~1 Fax : +82-42-934-5562

The World's best high-end Electromagnetic Field Therapy

Salus-Talent Pro

The World's 1st development 3 Tesla,
2Channel, Handheld Transducer
Magnetic Field Therapy.

Clinical Application

- * Acute or Chronic Back pain/Cervical pain, Sciatica
- * Osteoarthritis, rheumatoid Arthritis pain
- * Exercise/sensory nerve pain
- * Frozen Shoulder/neck pain/pelvic pain
- * Muscle Relaxation
- * Manipulation therapy by using hand held transducer

Large Transducer Application



Mini Transducer Application



CR-3000 (2003Year), Salus-Talent (2008Year).

Innovation Once Again


Salus-Talent Pro has smoother & stronger output for treating of deep sites by hand held Transducer.



Mechanism and Features



Principle of the Magnetic Stimulator

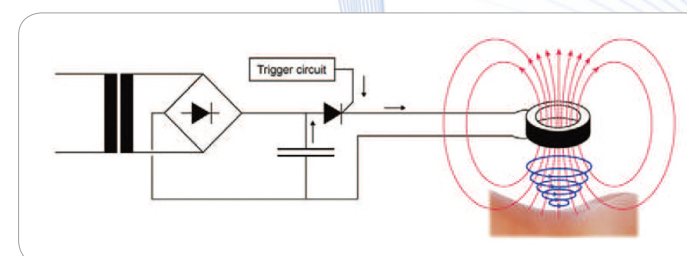

It generates strong magnetic field (applying eddy current deep inside) that can pass through clothes, tissues and bones which enables stimulation of a particular area deep inside the body without surgery. It almost does not cause pain during treatment and patients can get treatment with great ease.